



**T6: STEP BEYOND
EVERY EXPECTATION.**





Nearly 25 years ago, NuStep created the first product of its kind – the recumbent cross trainer. Now a mainstay in health, wellness, and fitness settings, the recumbent cross trainer has transformed hundreds of thousands of lives worldwide.

NuStep recumbent cross trainers offer an ideal workout solution for users who are fit, deconditioned, living with chronic conditions or recovering from a health event.

Today, NuStep continues its mission to transform lives with the launch of T6. Built on our proven T5 cross trainer platform, the T6 delivers the same cardio and strength workout users rely on, but with advanced features and programs to enhance their workout experience.


Features

Every detail was considered.

- > 15 levels of resistance, easily changed during your workout using the lower screen
- > Ergonomic hand grips rotate 40 degrees and sturdy arm handles, with easy-to-use clamshell release levers, adjust up to 18"
- > 12-degree seatback recline
- > Durable, oversized foot pedals have non-slip inserts and adjustable foot straps
- > Integrated media shelf
- > Compatible with T5 series accessories and adaptive equipment
- > Easy start up – system wakes up by touching the screen or stepping
- > Water bottle holder included
- > Optional grab ring and hip and torso belt for additional stability

A smart, sensible interface.

- > The upper screen provides an excellent interface for user engagement with an enhanced Pace Partner program, Quick Start and Charting Apps
- > The lower screen displays all workout metrics including time, speed (steps per minute), total distance, watts, calories, METs, resistance level and heartrate*
- > Admin mode allows for customization of the default settings
- > Heart rate monitoring through **Bluetooth®** wireless technology*
- > Audible feedback with volume control
- > 3-minute console shutdown with extended pause option for interval training
- > Seat position readout displayed on screen
- > USB port for data export to PDF
- > Compatible with NuStep App for iOS and Android™

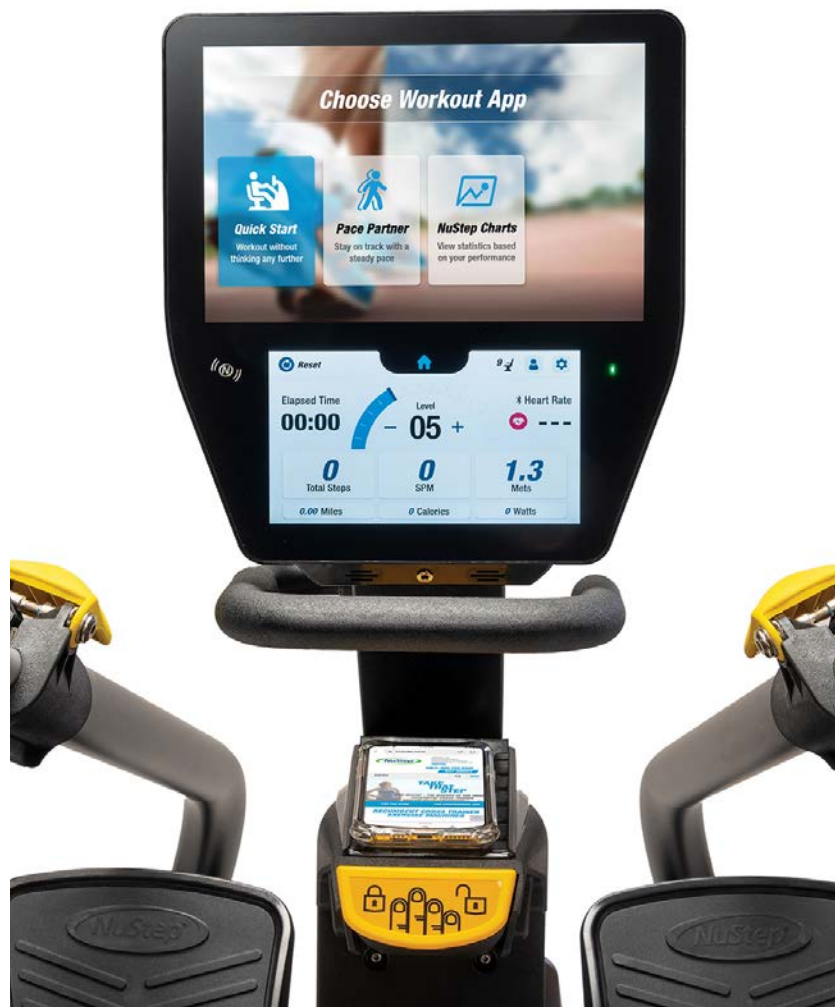

*When a Bluetooth heart rate monitor is connected.



Easy arm adjustments

40° of hand grip rotation

Dual-screen console with easy-to-read touchscreen display

Comfortable 360-degree swivel seat supports users up to 500 lbs

Low step-through design



StrideLock® feature locks handles and pedals



12-degree seatback recline offers maximum comfort



Adjustable foot straps to fit every user

Specs

Machine footprint: 73"L x 30"W x 52"H (185cm x 76cm x 117cm)

Weight: 285 lbs. (129 kg)

Standard seat width: 18" (46cm)

Wide seat width: 22" (56cm)

Fits users from 4'6" to 6'7" (137cm to 200cm) in height

Power: Medical grade 24V DC power supply included

Build your T6

Base Model Part No. Description

06001

T6_{PRO} with standard 18" width seat

06003

T6_{MAX} with 22" width seat

Optional Add-ons*

Part Number

Description

52341

Grab ring support

52421

Lap belt

59095

Torso belt

Adaptive Equipment

Part Number

Description

52000

Unilateral Leg Stabilizer Kit

Includes 1 leg stabilizer,

1 strap (26") and 2 storage clips

52065

Bilateral Leg Stabilizer Kit

Includes 2 leg stabilizers,

2 leg straps (26") and 2 storage clips

50680

Wellgrip® – left hand only

50681

Wellgrip – right hand only

50685

Wellgrip Set – left and right hand

Accessories

Part Number

Description

50780

Accessory Storage Bag

41815

Floor Mat

05300

NuStep Transporter

51877

Polar®H7 Heart Rate Sensor & Strap

*Must be ordered at time of T6 purchase. Not available as a retrofit option.

Why NuStep?

Quality:

We manufacture NuStep cross trainers at our headquarters in Ann Arbor, Michigan to ensure they meet strict quality standards.

Reputation:

Our superior products and longevity in the marketplace have earned us a reputation for excellence among industry professionals in healthcare, senior living and fitness.

Satisfaction:

We back every T6 cross trainer with a 30-day money back guarantee and a 10-year frame, 3-year parts and 1-year labor warranty.

Service:

Our customer service team – staffed by experienced product specialists – achieves an average 98% satisfaction rating among the customers they've assisted.

ORTHOCANADA